



**SMOBIO**

Small Bio, Smart Tool

[www.smobio.com](http://www.smobio.com)

## Product Information

### ExcelDye 6X DNA Loading Dye

**DL1000      5 ml × 2 (Orange, 6X)**

## Storage

4°C      ≥ 12 months

-20°C    ≥ 36 months

## Working Reagent Preparation

1:6 dilution in DNA electrophoresis sample.

(1µl loading dye for every 5 µl DNA sample)

## Description

SMOBiO ExcelDye 6X DNA Loading dye (Green) is pre-mixed buffer for tracking the DNA sample during the electrophoresis on agarose or polyacrylamide gels. It contains **Orange G** dye for tracking the DNA migration. The **Orange G** migrate at approximately 30 bp on a standard 2% TAE agarose gel (50 bp on 1% TAE agarose gel). The included glycerol can keep the DNA at the bottom of the well and the presence of EDTA chelates divalent metal ions to prevent the process of metal-dependent nuclease.

## Composition


0.15% Orange G  
10 mM Tris-HCl (pH 8.0)  
60% glycerol  
60 mM EDTA





## ExcelBand DNA Ladders

DNA Mass (ng/ 5µl) / Base Pairs (bp)

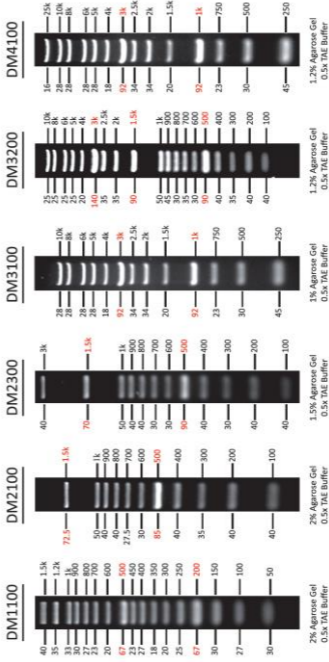



## FluoroBand DNA Ladders

